

A ROBUST FETAL ECG DETECTION METHOD FOR ABDOMINAL RECORDINGS

Biomedical Signal Processing Project

Niccolò Pizzocri - 75095A

INTRODUCTION


Monitoring the fetal cardiac activity can provide important information to obstetricians for the assessment of the **fetal well-being**.

Abdominal FECG provide a **non-invasive** diagnostic tool to do so.

However, the **fetal signal** has **small amplitude** and it's **overwhelmed by** a large number of **interference signals** such as baseline drift, power-line interference, motion artifacts and the mother's ECG itself

The paper proposes a non-blind method for FECG detection in abdominal recordings during pregnancy and labour.

METHOD

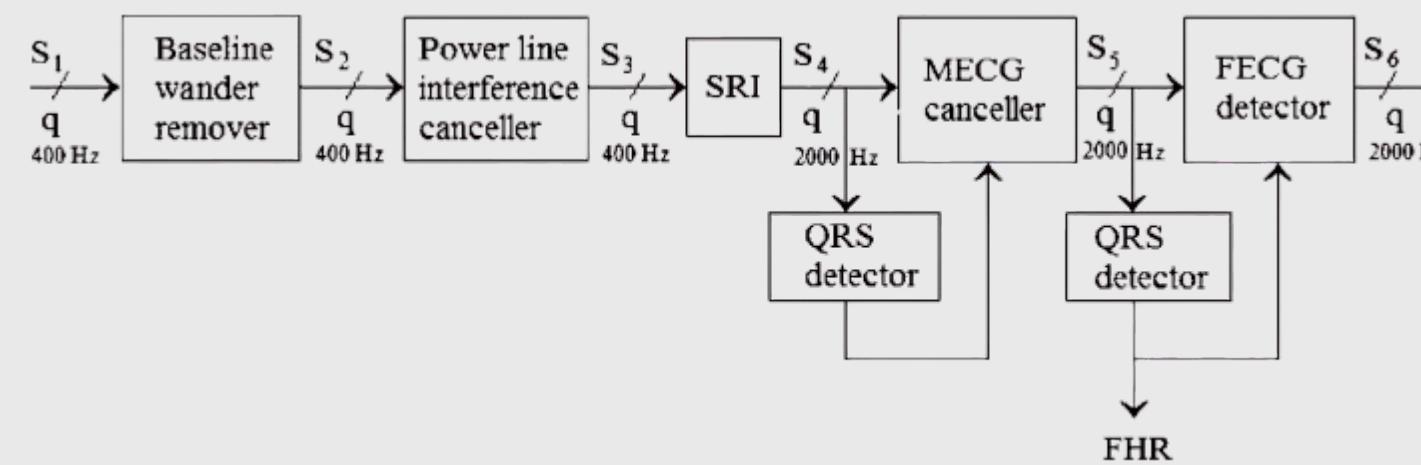


The method proposed detects the FECG by estimating and removing the inference signals step-by-step:

- 1.** Baseline wander remover
- 2.** Power-line interference canceller
- 3.** Sampling-rate increaser
- 4.** Mother QRS detector → MECG canceller
- 5.** Fetal QRS detector → FECG extractor

DATASET

2013 Computing in Cardiology Challenge's dataset:
One-minute of fetal ECG Recordings that include
four non-invasive abdominal signals
1000 Hz sampling rate

The dataset also contains the locations of each fetal QRS produced

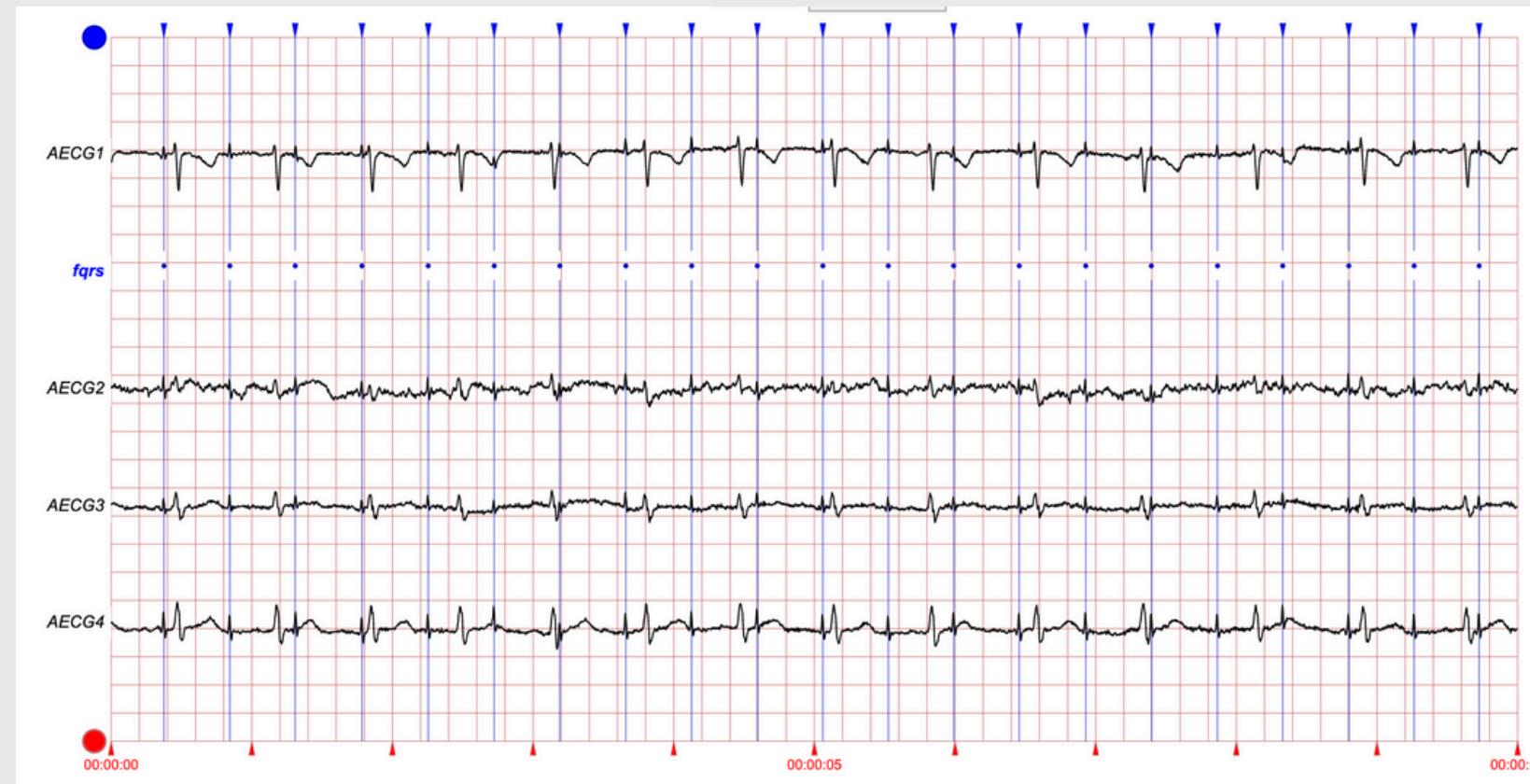

BASELINE WANDER REMOVER

Adaptive filtering and Polynomial or Cubic Spline Subtraction are rejected



SOLUTION → Linear-phase high pass FIR filter

It attenuates the large amount of baseline wander though it may cause a slight FECG distortion

Cut-off frequency: 3Hz

Filter Order: 1000 (+1)


Window: Hamming

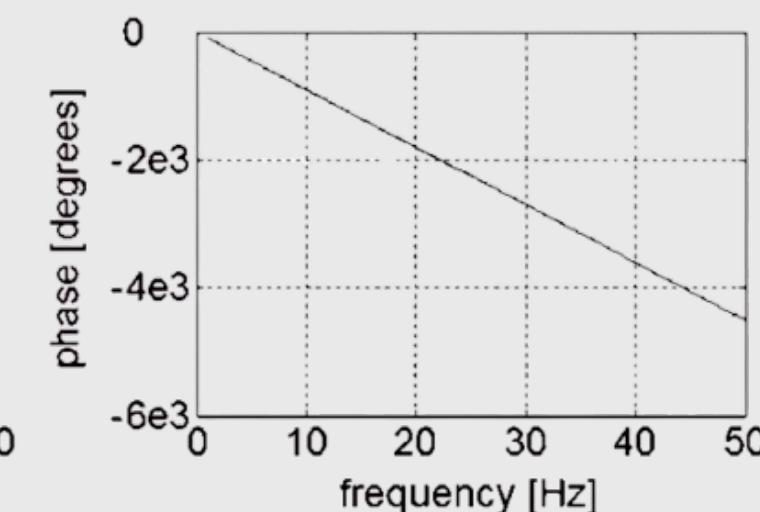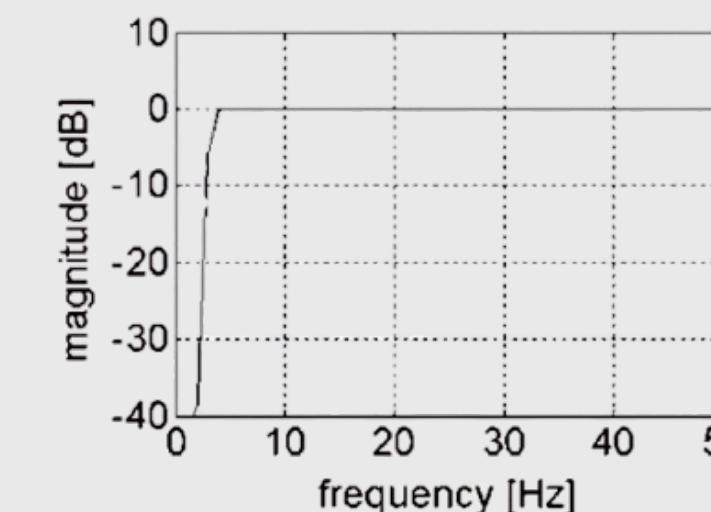

ADAPTIVE PLI CANCELLER


Blocking detection to protect QRS complexes:
Comb Filtering → Threshold Detection → Mask Generation

Adaptive Cancellation Loop:
Reference Generation → Error Calculation → Error Filtering
Reference Generation parameters Update using LMS

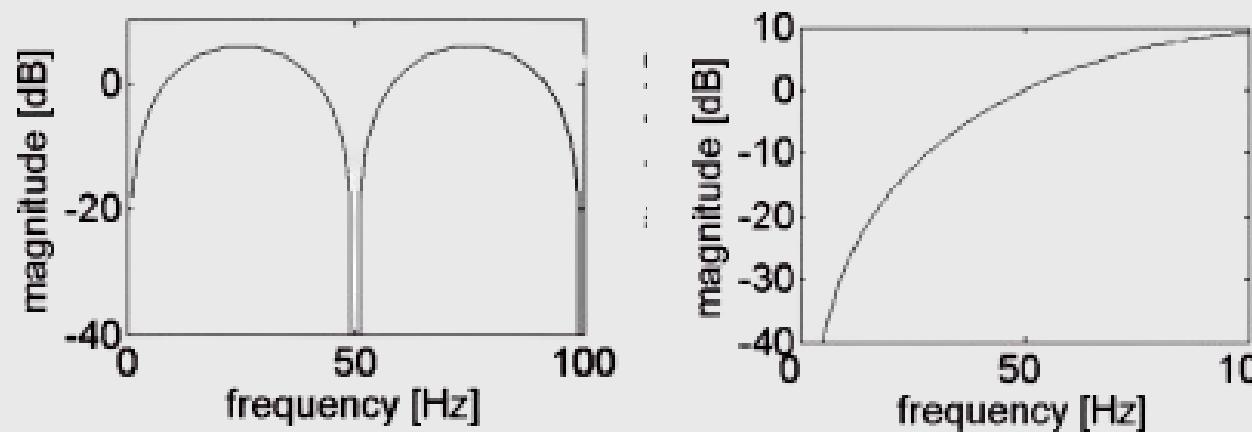

SAMPLING RATE INCREASER


Sampling rate increased to **2000 Hz** for removing MECG with high precision

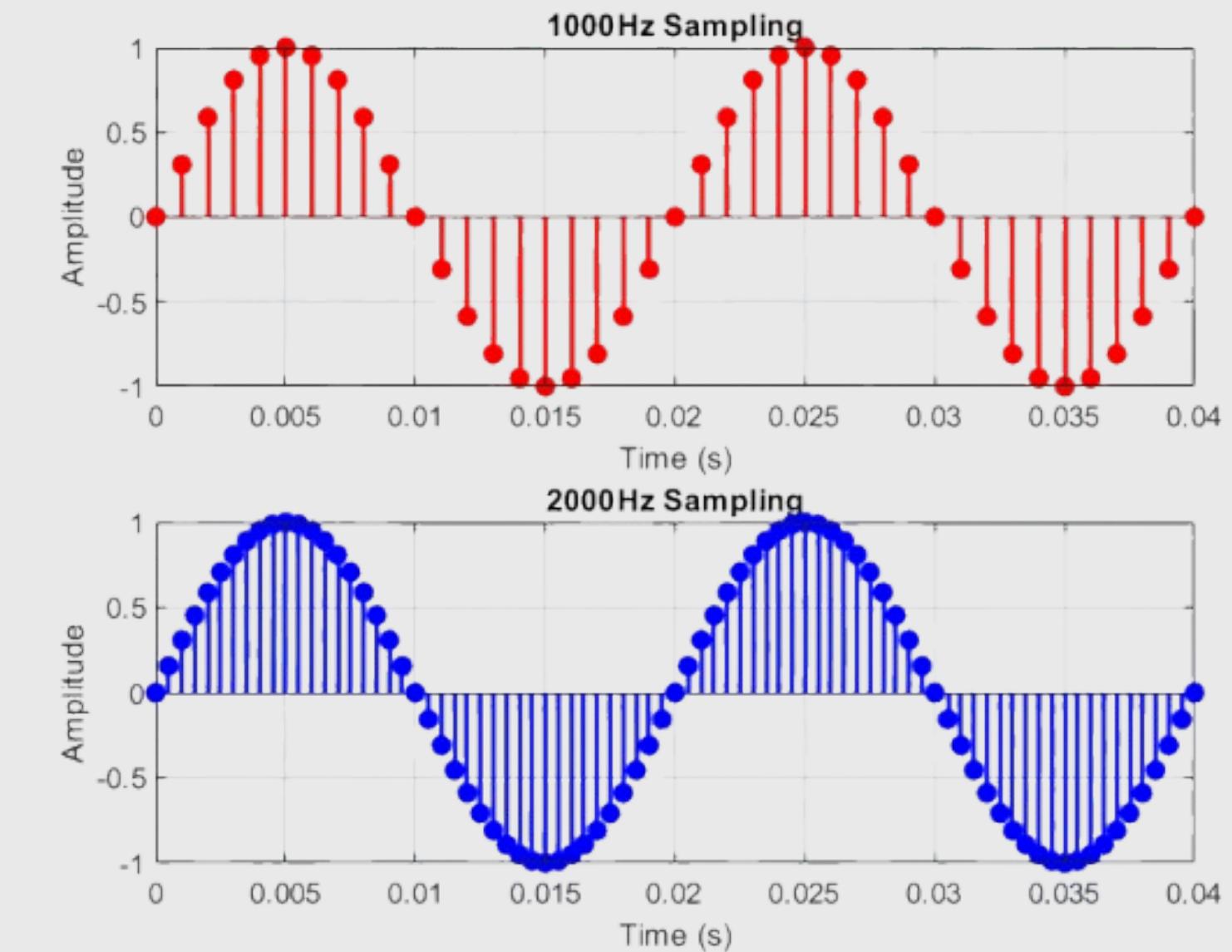

MECG CANCELLER

MQRS DETECTOR:

1. PCA
2. Detection with Hilber Envelope and Dynamic Thresholding
3. Perfect Template Overlap via Cross-Correlation

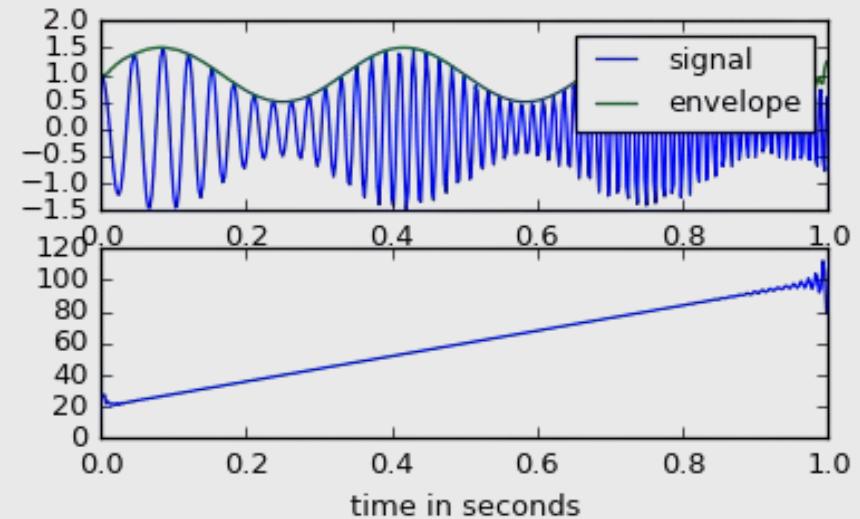

MECG CANCELLER:

1. Average Beat Template
2. Beat Segmented into P Wave, QRS complex and T wave
3. Adaptive Fitting with Ridge Regression

$$M = \begin{pmatrix} | & 0 & 0 \\ \underline{\mu}_P & 0 & 0 \\ | & 0 & 0 \\ 0 & | & 0 \\ 0 & \underline{\mu}_{QRS} & 0 \\ 0 & | & 0 \\ 0 & 0 & | \\ 0 & 0 & \underline{\mu}_T \\ 0 & 0 & | \end{pmatrix}. \quad \underline{a} = (M^T M)^{-1} M^T \underline{m}.$$

FECCG EXTRACTOR

FQRS DETECTOR

2 different approaches:


1. PCA-Hilbert
2. Energy-Based

ENERGY BASED DETECTION


- Compute Cubic Energy to enhance peaks
- Percentile Threshold

FETAL BEAT RECONSTRUCTION


- Synchronization
- Averaging of ALL beats

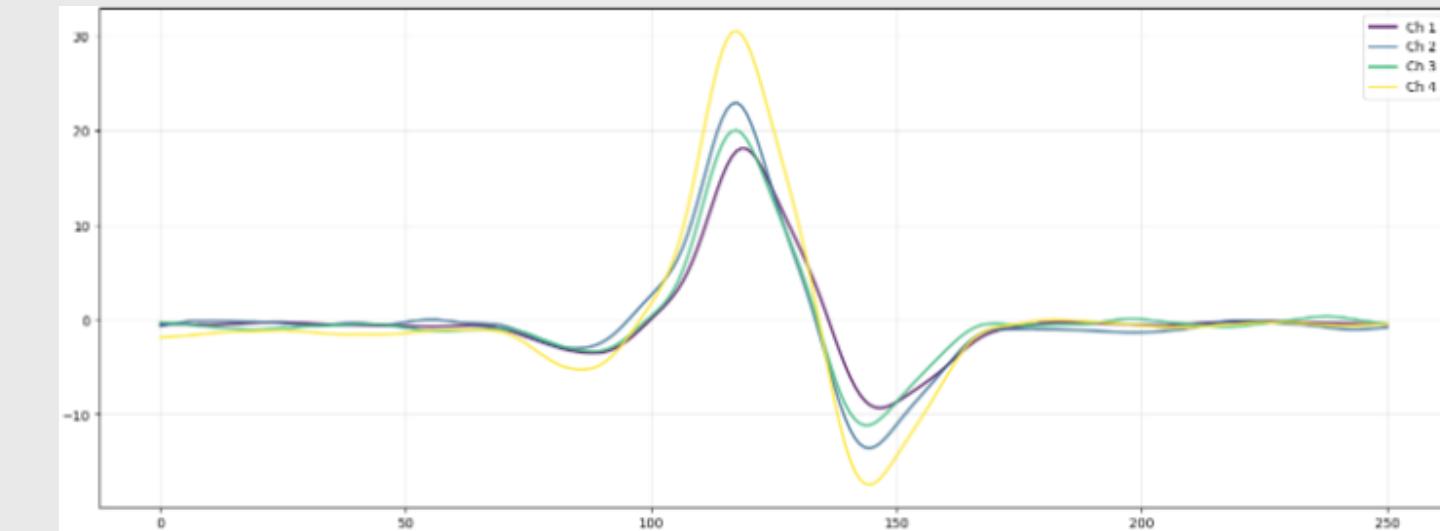

FECCG EXTRACTOR

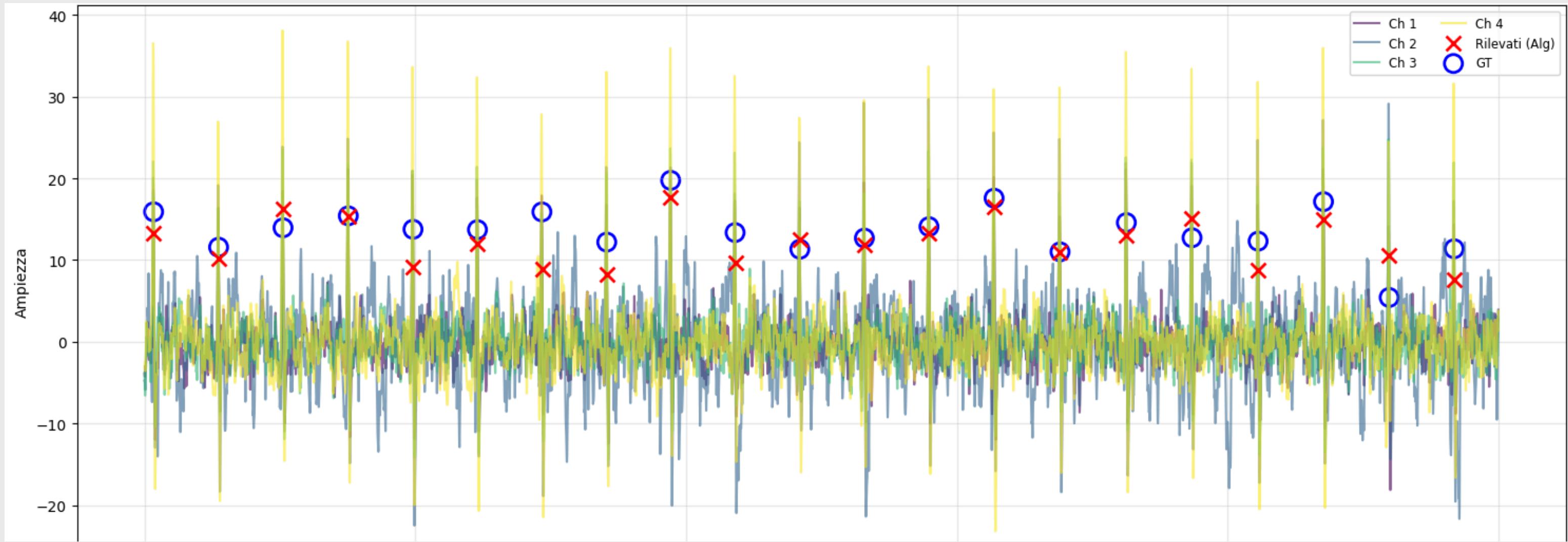

PIPELINE OUTPUT

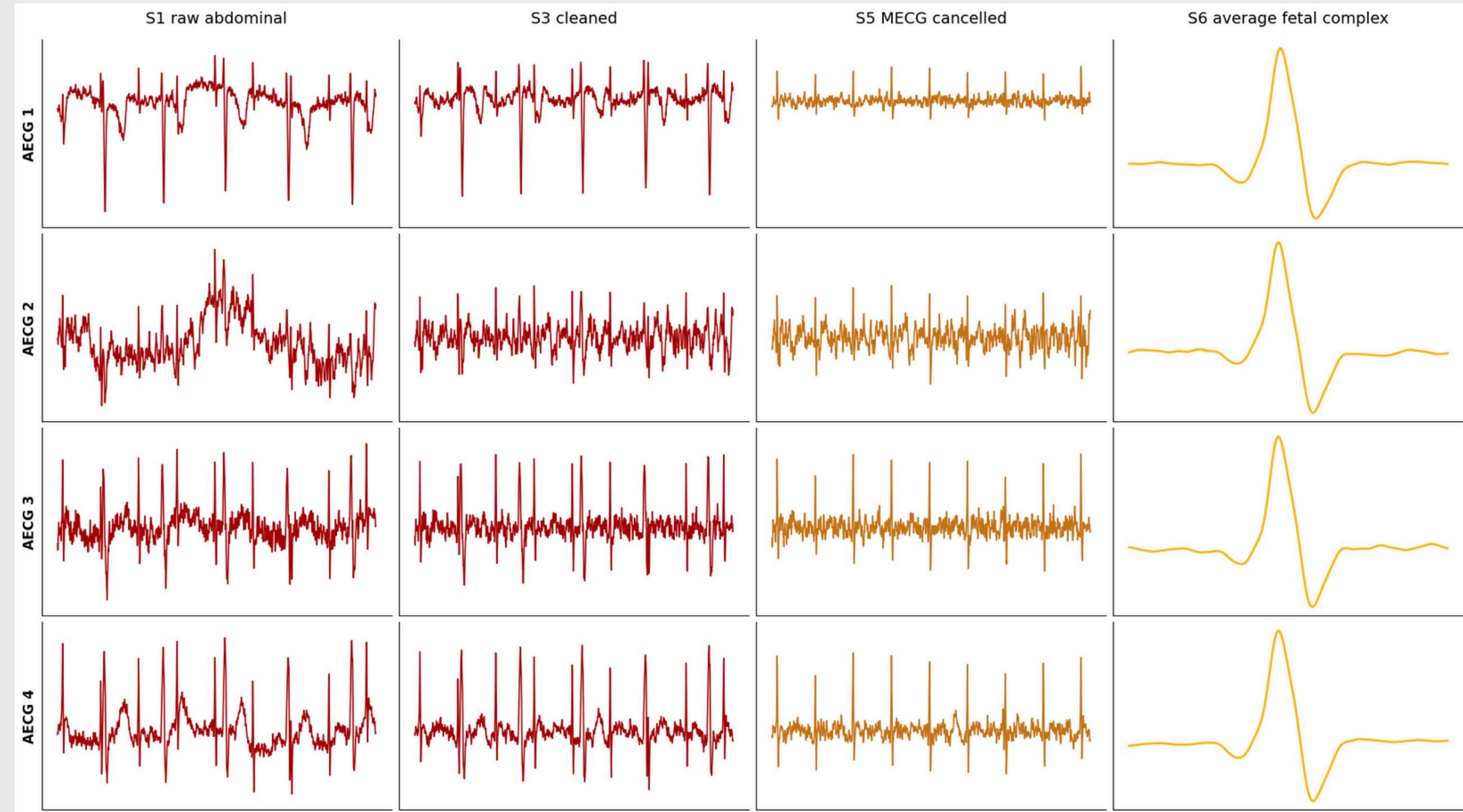

METRICS

$$Reliability = 1 - \frac{N_{outliers}}{N_{total}}$$

$$\text{SNR}_j = \frac{P_{F_j}}{P_{M_j}} = \frac{P_{S_{6,j}}}{P_{S_{5,j} - S_{6,j}}},$$

$$\text{SIR}_j = \frac{P_{F_j}}{P_{N_j}} = \frac{P_{S_{6,j}}}{P_{S_{4,j} - S_{5,j}}}.$$


$$\text{Accuracy} = \frac{TP + TN}{TP + TN + FP + FN}$$

$$\text{Sensitivity} = \frac{TP}{TP + FN}$$


$$\text{PPV (Precision)} = \frac{TP}{TP + FP}$$

$$\text{F1 Score} = \frac{2 \cdot TP}{2 \cdot TP + FP + FN}$$


RESULTS

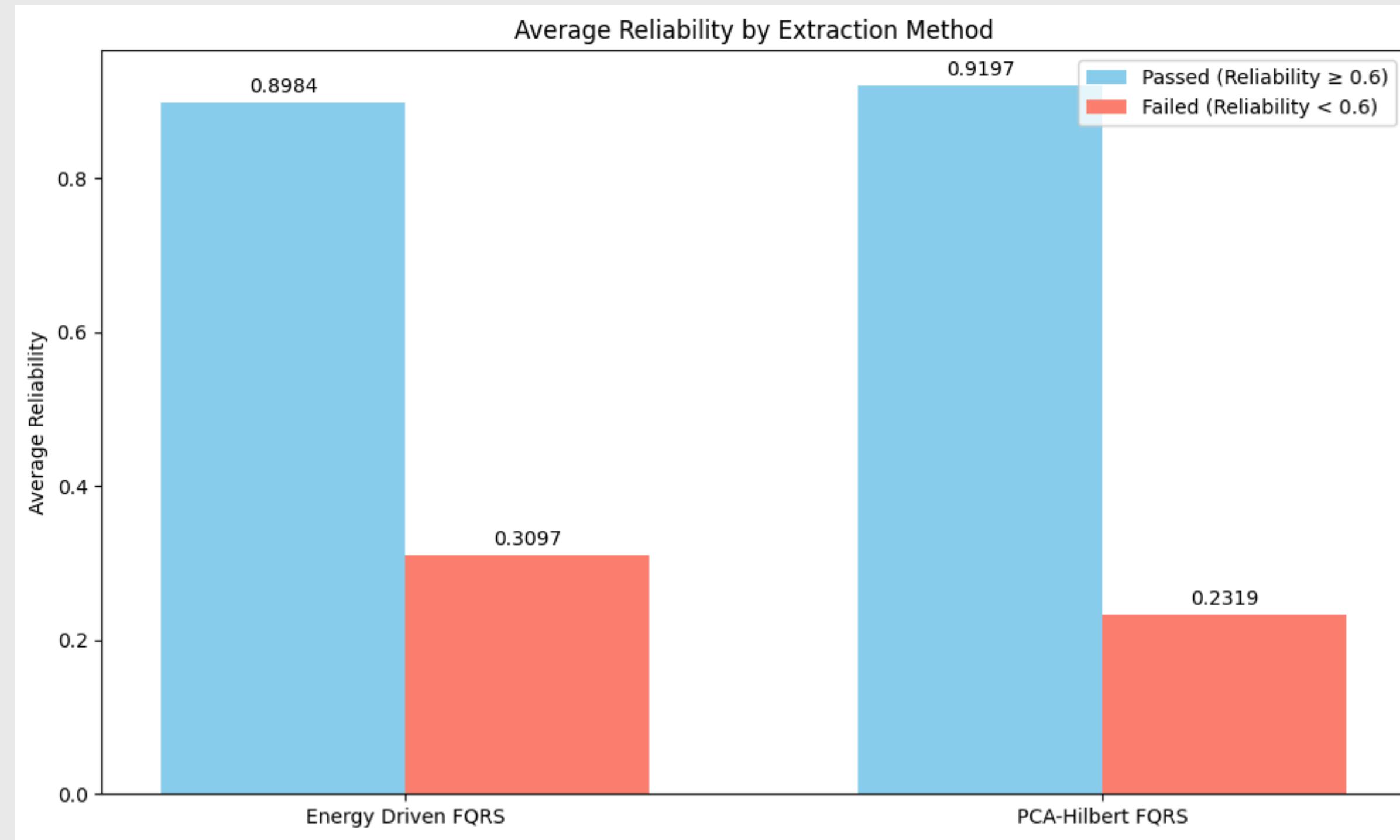

RESULTS

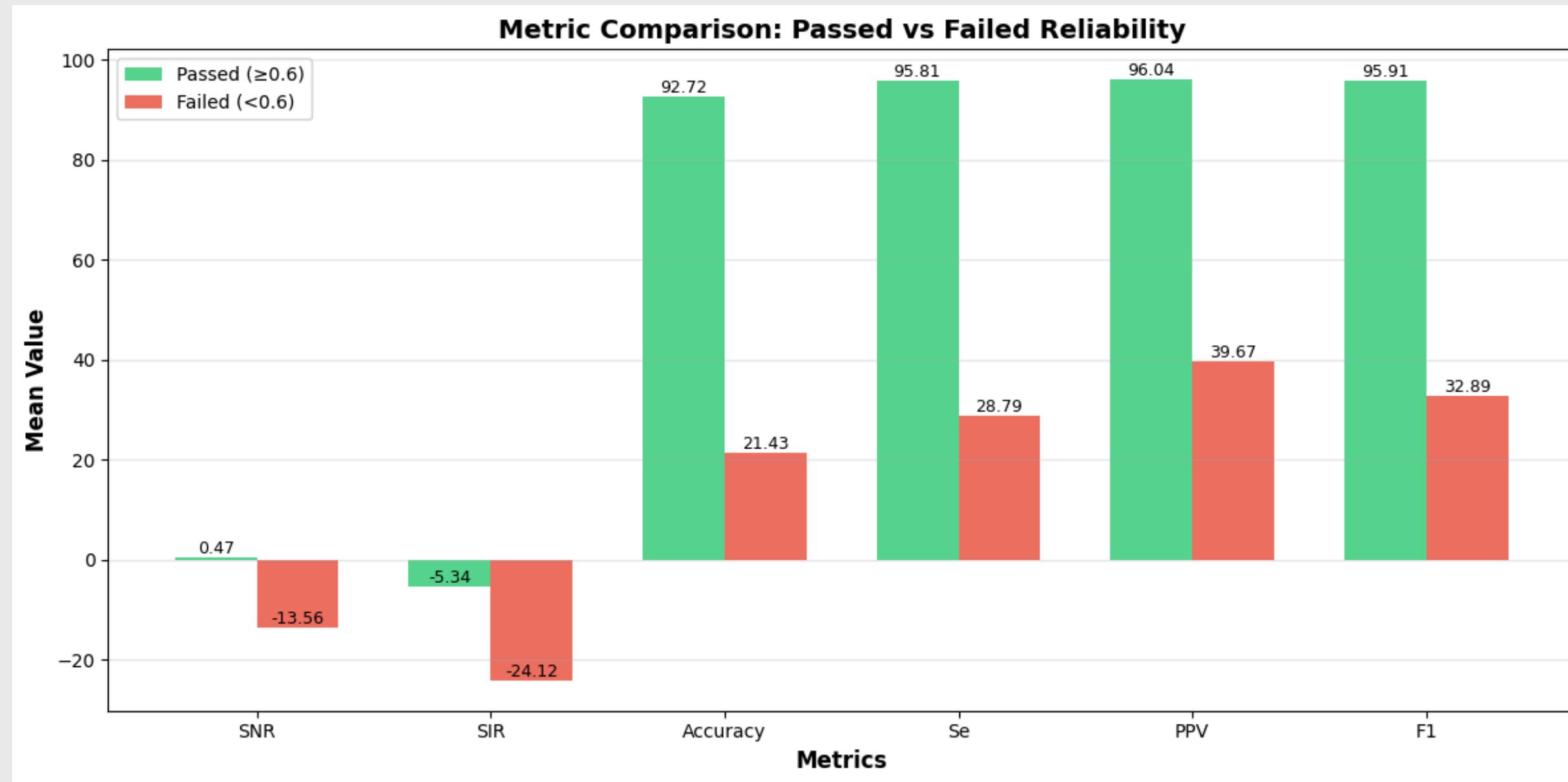


RESULTS

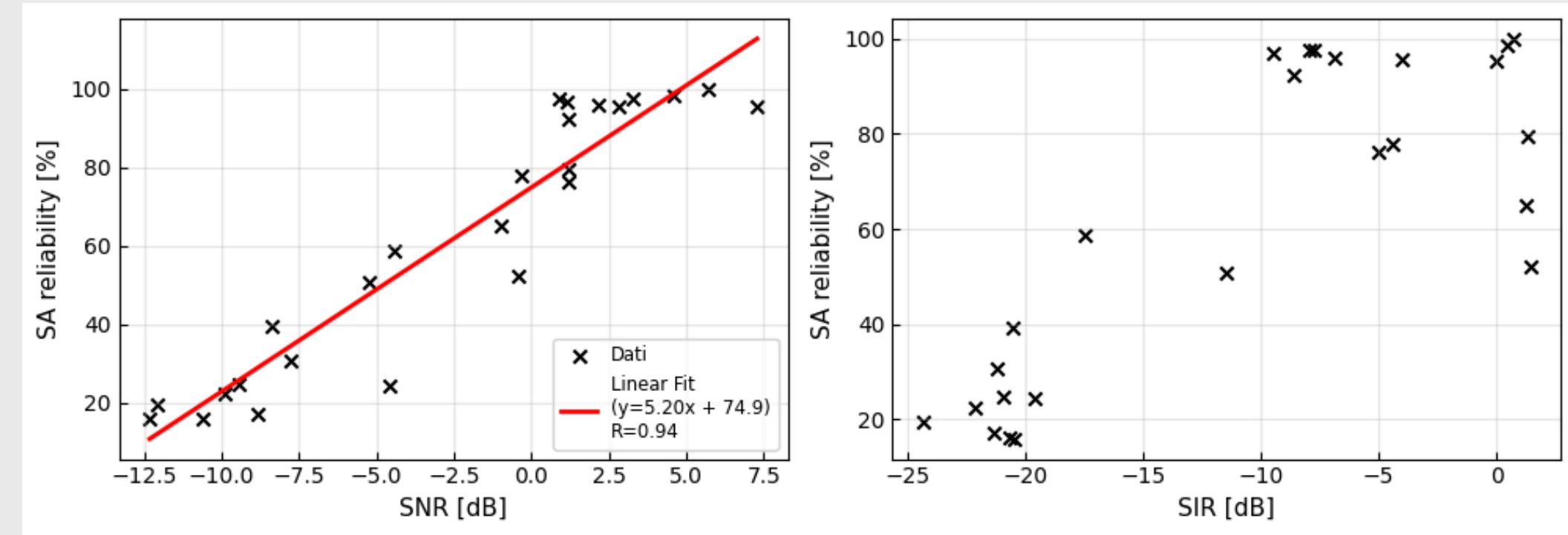


PEARSON CORRELATION MATRIX

	Reliability	SNR	SIR
Reliability	1	0.94	0.81
SNR	0.94	1	0.89
SIR	0.81	0.89	1

ICA COMPARISON

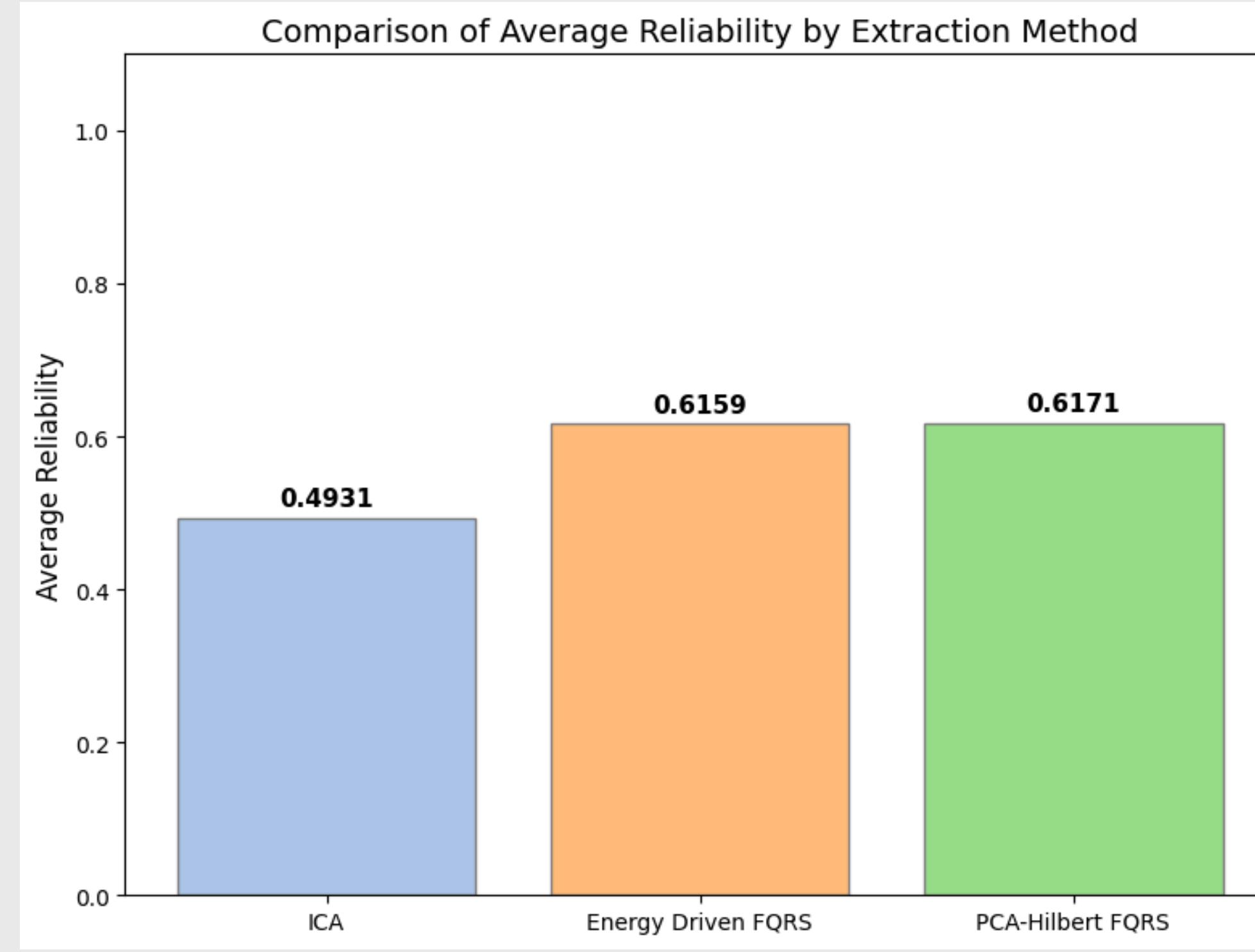

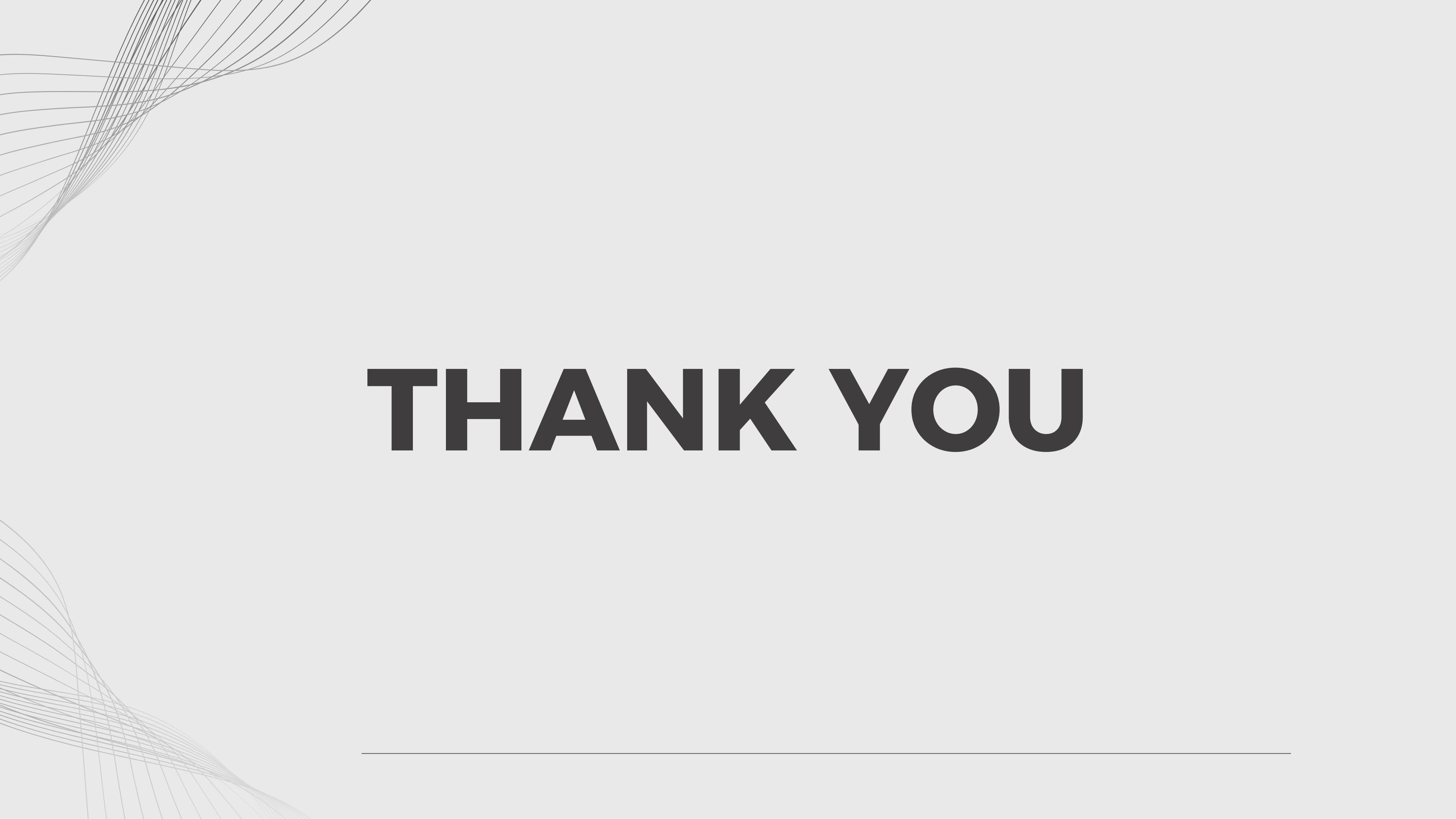

THANK YOU

